

PROOFREADING APPROVAL/RELEASED ARTWORK		Check as applicable
Printed Name		<input type="checkbox"/> Excluding dimensions
Signature	Date	<input type="checkbox"/> Excluding bar code dimensions
		<input type="checkbox"/> U.S. English
		<input type="checkbox"/> Other Language:
Printed Name		<input type="checkbox"/> Excluding dimensions
Signature	Date	<input type="checkbox"/> Excluding bar code dimensions
		<input type="checkbox"/> U.S. English
		<input type="checkbox"/> Other Language:



FLOSEAL Hemostatic Matrix, 5 mL

Instructions for Use

Caution: Federal Law (United States) restricts this device to sale by or on the order of a licensed healthcare practitioner.

DO NOT INJECT INTRAVASCULARLY.

FLOSEAL Hemostatic Matrix, also referred to as FLOSEAL Matrix, must not be injected into blood vessels.

Device Description:

The FLOSEAL Matrix kit consists of a bovine-derived Gelatin Matrix component, a human-derived Thrombin component, Applicator tips, and several mixing accessories. The mixing accessories include a syringe with an integral female Luer connector, a Thrombin bowl, a pre-filled Sodium Chloride Solution syringe and a vial adapter for needle-free reconstitution. The accessories are included to facilitate the reconstitution and mixing of the Thrombin into the Gelatin Matrix. Applicator tips are included to facilitate the delivery of FLOSEAL Matrix to the site to be treated. (For specific package contents, see Table in "How Supplied" section.)

The Gelatin Matrix consists of cross-linked gelatin granules and is provided sterile and non-pyrogenic in a standard disposable syringe. The Thrombin (Human) is a sterile, non-pyrogenic, freeze-dried, vapor-heated and solvent/detergent treated powder preparation made from pooled human plasma. The Sodium Chloride Solution is a sterile, non-pyrogenic solution. After reconstitution of the lyophilized Thrombin with Sodium Chloride Solution, the resulting Thrombin solution contains 500 units/mL Thrombin (Human).

*The potency expressed in units is determined using a clotting assay against an internal reference standard for potency that has been calibrated against the World Health Organization (WHO) Second International Standard for Thrombin, 01/580. Therefore, a unit used herein is equivalent to an International Unit.

Thrombin (Human) is prepared from pooled human plasma through a series of separation and filtration steps followed by incubation of the solution with Calcium Chloride in order to activate prothrombin to Thrombin. The Thrombin (Human) solution subsequently undergoes ultra/diafiltration, vapor heat treatment, solvent/detergent treatment, sterile filtration and freeze-drying.

Thrombin (Human) is made from pooled human plasma obtained from U.S. licensed plasma collection centers. The two-step vapor-heated and solvent/detergent treatment used in its manufacture has been shown to be capable of significant viral reduction. However, no procedure has been shown to be completely effective in removing viral infectivity from derivatives of human plasma (see "Warnings").

The manufacturing procedure for FLOSEAL Matrix includes processing steps designed to reduce the risk of viral transmission. Several steps are included in the manufacture of the Gelatin Matrix component that reduce the risk of viral transmission. The virus reduction factors (expressed as log₁₀) for the manufacture of the Gelatin Matrix component are provided in the table below.

Reduction Factors for Virus Removal and/or Inactivation During the Manufacture of Gelatin Matrix		
Manufacturing Step	Virus Reduction Factor of Virus Tested	
	BVDV	PPV
Base Treatment (NaOH)	>5.4	4.0
Chemical Cross-linking	>5.0	1.1
Heat Treatment	>6.5	1.9

BVDV: Bovine viral diarrhea virus, a model for Hepatitis C virus; **PPV:** Porcine parvovirus, a model for non-enveloped viruses, among those, Hepatitis A virus.

A two-step vapor-heated and solvent/detergent viral inactivation treatment process is included in the manufacture of Thrombin. The virus reduction factors (expressed as log₁₀) for Thrombin are provided in the table below.

Reduction Factors for Virus Removal and/or Inactivation Thrombin Component						
Manufacturing Step	Mean Reduction Factors (log ₁₀) of Virus Tested					
	HIV-1	HAV	BVDV	PRV	MMV	B19V
Thrombin precursor mass capture	3.2	1.5	1.8	2.5	1.2	1.7
Vapor Heat Treatment	>5.5	>4.9	>5.3	>6.7	1.0	>4
Solvent/Detergent Treatment	>5.3	n.d.	>5.5	>6.4	n.d.	n.d.
Ion Exchange Chromatography	n.d.	n.d.	n.d.	n.d.	3.6	n.d.
Overall Reduction Factor (ORF)	>14.0	>6.4	>12.6	>15.6	5.8	>5.7

n. d. = not determined

HIV-1: Human immunodeficiency virus 1; **HAV:** Hepatitis A virus; **BVDV:** Bovine viral diarrhea virus, a model for Hepatitis C virus; **PRV:** Pseudorabies virus, a model for enveloped DNA viruses, among those Hepatitis B virus; **MMV:** Mice minute virus, a model for Human Parvovirus B19; **B19V:** Human Parvovirus B19.

FLOSEAL Matrix is the combination of the Gelatin Matrix component and the reconstituted Thrombin (Human) component. Thrombin must be added to the Gelatin Matrix component prior to use. FLOSEAL Matrix is biocompatible and resorbed within 6 to 8 weeks, consistent with normal wound healing. FLOSEAL Matrix is intended only for topical administration.

Indications:

FLOSEAL Matrix is indicated in surgical procedures (other than in ophthalmic) as an adjunct to hemostasis when control of bleeding, by ligature or conventional procedures is ineffective or impractical.

Contraindications:

- Do not inject or compress FLOSEAL Matrix into blood vessels. Do not apply FLOSEAL Matrix in the absence of active blood flow, e.g., while the vessel is clamped or bypassed, as extensive intravascular clotting and even death may result.
- To avoid a risk of allergic-anaphylactoid reaction and/or thromboembolic events, which may be life-threatening, do not inject FLOSEAL Matrix into a vessel or tissue.

- Do not use FLOSEAL Matrix in patients with known allergies to materials of bovine origin.
- Do not use FLOSEAL Matrix in the closure of skin incisions because it may interfere with the healing of the skin edges due to mechanical interposition of gelatin.

Warnings:

- FLOSEAL Matrix contains Thrombin made from human plasma. Products made from human plasma may contain infectious agents, such as viruses, that can cause disease. The risk that such products will transmit an infectious agent has been reduced by screening plasma donors for prior exposure to certain viruses, by testing for the presence of certain current virus infections, and by inactivating and removing certain viruses. Despite these measures, such products can still potentially transmit disease. Because this product is made from human blood, it may carry a risk of transmitting infectious agents, e.g., viruses, and theoretically, the Creutzfeldt-Jakob disease (CJD) agent. All infections thought by a physician possibly to have been transmitted by this product should be reported by the physician or other healthcare provider to Baxter Healthcare Corporation. The physician should discuss the risks and benefits of this product with the patient.
- FLOSEAL Matrix is not intended as a substitute for meticulous surgical technique and the proper application of ligatures or other conventional procedures for hemostasis.
- FLOSEAL Matrix is not intended to be used as a prophylactic hemostatic agent.
- Excess FLOSEAL Matrix (material not incorporated in the hemostatic clot) should always be removed by gentle irrigation from the site of application. Meticulous irrigation is required when used in, around, or in proximity to foramina in bone, areas of bony confine, the spinal cord, the brain and/or cranial nerves.
- FLOSEAL Matrix swells by approximately 10-20% after product is applied and surgeons should consider its potential effect on the surrounding anatomic areas. Maximum swell volume is achieved within about 10 minutes.
- FLOSEAL Matrix should not be used in the presence of infection. FLOSEAL Matrix should be used with caution in contaminated areas of the body. If signs of infection or abscess develop where FLOSEAL Matrix has been applied, re-operation may be necessary in order to remove the infected material and allow drainage.
- The safety and effectiveness of FLOSEAL Matrix for use in ophthalmic procedures has not been established.
- FLOSEAL Matrix should not be used for controlling post-partum bleeding or menorrhagia.
- The safety and effectiveness of FLOSEAL Matrix has not been established in children under 2 years of age and pregnant women.
- Do not use air to remove residual FLOSEAL Matrix from Applicator tip.

Precautions:

General

- For single use only. Do not re-sterilize.
- As with other hemostatic agents, circumstances that result in a negative peripheral venous pressure (e.g. patient positioning) may draw material into the vascular system, potentially resulting in life-threatening thromboembolic events.
- Since the Thrombin solution can be denatured by contact with solutions containing alcohol, iodine, or heavy metal ions, FLOSEAL Matrix should not be applied before the application site is cleaned to remove any antiseptics that may contain such substances.
- When placed into cavities or closed tissue spaces, gentle approximation is advised.
- As with other hemostatic agents, do not aspirate FLOSEAL Matrix into extracorporeal cardiopulmonary bypass circuits or autologous blood salvage circuits. It has been demonstrated that fragments of collagen based hemostatic agents may pass through 40µm transfusion filters of blood scavenging systems.
- Do not use FLOSEAL Matrix on bone surfaces where adhesives, such as methylmethacrylate or other acrylic adhesives, will be required to attach a prosthetic device. Microfibrillar collagen has been reported to reduce the strength of methylmethacrylate adhesives used to attach prosthetic devices to bone surfaces.
- FLOSEAL Matrix should not be used for the primary treatment of coagulation disorders.
- The safety and effectiveness of the combined use of FLOSEAL Matrix with antibiotic solutions or powders has not been established.
- The safety and effectiveness for use in neurosurgical and urological procedures has not been established through randomized clinical studies.
- In urological procedures, FLOSEAL Matrix should not be left in the renal pelvis or ureters to eliminate the potential foci for calculus formation.
- The Applicator tips should not be cut.

Information for Patients

- Some viruses, such as human parvovirus B19, are particularly difficult to remove or inactivate at this time. Human parvovirus B19 most seriously affects pregnant women, or immune-compromised individuals. Symptoms of human parvovirus B19 infection include fever, drowsiness, chills, and runny nose followed about two weeks later by a rash, and joint pain. Patients should be encouraged to consult their physician if such symptoms appear.

Carcinogenesis, Mutagenesis, Impairment of Fertility

- Long-term animal studies to evaluate the carcinogenic potential of FLOSEAL Matrix or studies to determine the effect of FLOSEAL Matrix on fertility have not been performed.

Use in Pregnancy

- It is not known whether FLOSEAL Matrix can cause fetal harm when administered to a pregnant woman or can affect reproduction capacity. FLOSEAL Matrix should be administered to a pregnant woman only if clearly needed.

FLOSEAL Matrix Adverse Events:

In a randomized prospective, concurrently controlled clinical trial using a formulation of FLOSEAL Matrix containing bovine Thrombin, a total of 309 patients received FLOSEAL Matrix or the Control (Gelatin Sponge + Thrombin). The most common adverse events recorded during and after the application of the hemostatic agents were anemia, atrial fibrillation, infection, and hemorrhage. The following is a complete list of adverse events reported in greater than 1% of patients that were observed in the pivotal clinical trial for the FLOSEAL Matrix group. The corresponding adverse events for the Control group are listed for comparison. None of the adverse events that occurred were judged by the surgeon to be "Probably Related" to the use of FLOSEAL Matrix.

Adverse Events Reported in Greater than 1% of Patients in the FLOSEAL Matrix Clinical Trial		
Adverse Event	FLOSEAL Matrix	Control (Gelatin Sponge + Thrombin)
Anemia	12 (8%)	7 (4%)
Fibrillation Atrial	10 (6%)	8 (5%)
Infection	10 (6%)	11(7%)
Hemorrhage	6 (4%)	6 (4%)
Pneumonia	6 (4%)	2 (1%)
Urinary Tract Infection	6 (4%)	3 (2%)
Rash	5 (3%)	3 (2%)
Edema	5 (3%)	1 (<1%)
Hypotension	4 (3%)	2 (1%)
Respiratory Distress	4 (3%)	3 (2%)
Confusion	4 (3%)	0 (0%)
Dural Tear	4 (3%)	4 (3%)
Fibrillation Ventricular	4 (3%)	3 (2%)
Arrhythmia	4 (3%)	0 (0%)
Heart Failure Right	3 (2%)	2 (1%)
Thrombosis Arterial	3 (2%)	8 (5%)
Fever	3 (2%)	2 (1%)
Atelectasis	3 (2%)	1 (<1%)
Pleural Effusion	3 (2%)	5 (3%)

Counts reflect number of patients in each treatment group reporting one or more adverse events that map to a Modified COSTART 5th edition body system. At each level of summarization (Adverse Event), patients are only counted once.

Other adverse events observed in 1% or less of the FLOSEAL Matrix clinical trial patients were myocardial infarction, cellulitis, pneumothorax, pain, cerebrovascular accident, hallucination, paresthesia, bradycardia, abscess, diarrhea, urinary retention, dehiscence, skin ulcer, transfusion reaction, dyspnea, heart arrest, lung edema, back pain, ventricular tachycardia, neuropathy, acute kidney failure, kidney tubule necrosis, gastritis, nausea, and vomiting, skin rash, hyperglycemia, and heel ulcer.

The following adverse events, all rated "mild", were deemed by the surgeon to be "Possibly Related" to the use of FLOSEAL Matrix: anemia (2 patients, 1%), mild post-operative bleeding (1 patient, <1%), and local inflammation (1 patient, <1%). No other adverse events were deemed by the surgeon to be related to the use of FLOSEAL Matrix.

Allergic reactions may be encountered in people who are sensitive to bovine materials.

Gelatin-Based Hemostatic Agents Adverse Events:

- Gelatin-based hemostatic agents may serve as a nidus for infection and abscess formation and have been reported to potentiate bacterial growth.
- Excess gelatin-based hemostatic agents should always be removed by gentle irrigation from the site of application. Removal of excess is done to avoid excessive inflammatory reaction, adhesion and/or granuloma formation.
- Giant cell granulomas have been observed at implant sites when used in the brain.
- Compression of the brain and spinal cord resulting from the accumulation of sterile fluid has been observed.
- Multiple neurologic events were reported when absorbable gelatin-based hemostatic agents were used in laminectomy operations, including cauda equina syndrome, spinal stenosis, meningitis, arachnoiditis, headaches, paresthesias, pain, bladder and bowel dysfunction, and impotence.
- The use of absorbable gelatin-based hemostatic agents during the repair of dural defects associated with laminectomy and craniotomy operations has been associated with fever, infection, leg paresthesias, neck and back pain, bladder and bowel incontinence, cauda equina syndrome, neurogenic bladder, impotence, and paresis.
- The use of absorbable gelatin-based hemostatic agents has been associated with paralysis, due to device migration into foramina in the bone around the spinal cord, and blindness, due to device migration in the orbit of the eye, during lobectomy, laminectomy and repair of a frontal skull fracture and lacerated lobe.
- Foreign body reactions, "encapsulation" of fluid, and hematoma have been observed at implant sites.
- Excessive fibrosis and prolonged fixation of a tendon have been reported when absorbable gelatin-based sponges were used in severed tendon repair.
- Toxic shock syndrome was reported in association with the use of absorbable gelatin-based hemostats in nasal surgery.
- Fever, failure of absorption, and hearing loss have been observed when absorbable hemostatic agents were used during tympanoplasty.

Human Thrombin Adverse Events:

As with any other plasma derivatives, anaphylactoid or anaphylactic reactions may occur in rare cases. No adverse events of this type were reported during the course of clinical trials using a different product containing the same human Thrombin component. Mild reactions can be managed with antihistamines; severe hypotensive reactions require immediate intervention using current principles of shock therapy.

Equivalence of Bovine and Human Thrombin:

The performance of FLOSEAL Matrix containing human Thrombin was compared to that of original FLOSEAL Matrix (containing bovine Thrombin) in a bleeding liver square model in pigs.

Blood flow rates for the lesions created in the pig liver model were recorded at specific time points and statistically analyzed by the method of Blackwelder and Chang modified for continuous variables. This analysis demonstrates that the performance of FLOSEAL Matrix (human Thrombin) is equivalent to the performance of FLOSEAL Matrix (bovine Thrombin) with a p-value of < 0.001 at each of the time intervals.

In addition, each lesion was subjectively scored for bleeding at each time point. These data points were analyzed using the method of Blackwelder and Chang for proportions. The results for all lesions in all animals showed FLOSEAL Matrix (human Thrombin) and FLOSEAL Matrix (bovine Thrombin) were equivalent for each of the time intervals with a p-value of 0.015.

Clinical Studies:

Study Design and Objectives: A prospective, randomized, controlled, multi-center, multi-specialty study was conducted using a formulation of FLOSEAL Matrix containing bovine Thrombin. Three hundred and nine (309) patients were enrolled at 10 centers. The objective of the study was to evaluate the safety and effectiveness of FLOSEAL Matrix, compared to a commercially available control hemostat, Absorbable Gelatin Sponge, USP ("Gelatin Sponge") + Thrombin, in controlling intraoperative bleeding. This study was designed to show that the FLOSEAL Matrix success rate was equivalent to the success rate for the Control.

Patients undergoing surgery in cardiac, vascular or spinal/orthopedic surgical specialties were included.

Patients were randomized only after it was determined that the bleeding could not be controlled using conventional approaches (e.g. direct pressure, sutures and/or cautery) because of their ineffectiveness or impracticality. Success at achieving hemostasis was defined as cessation of bleeding within 10 minutes following application of the agent. The primary endpoint was hemostasis success for the first treated bleeding site. A secondary endpoint was time to hemostasis for the first treated bleeding site. Although multiple bleeding sites in the same patient were treated, only the first treated bleeding site was used to determine primary effectiveness, as this was the only site that was truly randomized.

Clinical Study Results:

Primary Endpoint: The primary endpoint, cessation of bleeding within 10 minutes of the first lesion, achieved a success rate of 96% in the FLOSEAL Matrix group and 77% in the Control group. FLOSEAL Matrix and Control were shown to be equivalent using the Blackwelder and Chang test, using a Δ (clinically significant difference) of 0.15 (p<0.0001). The difference between FLOSEAL Matrix and Control was also shown to be statistically significant using the Cochran-Mantel-Haenszel test (p<0.001). Primary endpoint data were stratified for individual surgical specialties, and the results are summarized in the table below:

Hemostasis Within 10 Minutes – First Lesion Only (Intent-to-Treat Patients)		
Patient Category	FLOSEAL Matrix	Control
All Patients	96% (149/156)	77% (118/153)
Cardiac	94% (45/48)	60% (27/45)
Vascular	93% (40/43)	76% (35/46)
Spinal/Orthopedic	98% (64/65)	90% (56/62)

In the cardiac cohort, 88 of the 93 patients (95%) underwent surgery with extracorporeal cardiopulmonary bypass. FLOSEAL Matrix was used for hemostasis prior to heparin reversal by the administration of protamine sulfate in 19 of 46 patients. Protamine sulfate reverses the anticoagulative effects of heparin. Results for hemostasis at 10 minutes for the heparinized patients in both the FLOSEAL Matrix and Control groups, before and after protamine sulfate reversal of heparin, are shown in the table below:

Hemostasis Success at 10 Minutes Before and After Protamine Administration (Cardiac Patients Only)		
Group	Before Protamine	After Protamine
FLOSEAL Matrix	89% (17/19)	96% (26/27)
Control	36% (5/14)	75% (21/28)

The success rate for FLOSEAL Matrix did not appear to be affected by whether or not the patient had received protamine sulfate administration. This was demonstrated by the fact that the success rate for FLOSEAL Matrix before protamine sulfate administration was similar to the success rate after protamine sulfate administration whereas the Control hemostat success rate was clearly lower before protamine sulfate reversal of heparin was administered.

Secondary Endpoint: A secondary endpoint was time to hemostasis for the first treated bleeding site. The data for time to hemostasis are summarized in the table below.

Cumulative Percent of Patients with Complete Hemostasis First Lesion (Protocol Valid Patients*)		
Time Interval	FLOSEAL Matrix	Control
0 – 1 minute	41% (62/153)	21% (32/150)
1 – 2 minutes	69% (106/153)	32% (48/150)
2 – 3 minutes	85% (130/153)	48% (72/150)
3 – 6 minutes	93% (143/153)	68%(102/150)
6 – 10 minutes	97% (149/153)	77% (115/150)

**Six (6) patients, 3 in the FLOSEAL Matrix group and 3 in the Control group, were excluded because of protocol deviations in measuring hemostasis for the first treated bleeding site.*

When the data were stratified by surgical specialty, the median times to hemostasis were shorter for the FLOSEAL Matrix group than for the Control group in all specialties. The median times are summarized in the table below.

Time to Hemostasis First Lesion Only (Protocol Valid Lesions)		
Median Time to Hemostasis in minutes (95% Confidence Interval*)		
Patient Category	FLOSEAL Matrix	Control
All Patients	2.0 (1.5, 2.5)	6.0 (5.5, 6.0)
Cardiac	2.8 (2.0, 4.0)	8.0 (6.0, 8.5)
Vascular	2.5 (2.0, 4.0)	6.5 (4.5, 8.0)
Spinal/Orthopedic	1.5 (1.0, 1.5)	3.0 (2.0, 4.5)

*Confidence interval using a Bonferroni correction.

PROOFREADING APPROVAL/RELEASED ARTWORK		Check as applicable
Printed Name		<input type="checkbox"/> Excluding dimensions <input type="checkbox"/> Excluding bar code dimensions <input type="checkbox"/> U.S. English <input type="checkbox"/> Other Language:
Signature	Date	
Printed Name		<input type="checkbox"/> Excluding dimensions <input type="checkbox"/> Excluding bar code dimensions <input type="checkbox"/> U.S. English <input type="checkbox"/> Other Language:
Signature	Date	

Use of FLOSEAL Matrix as a Hemostatic Agent for Nasal/ Sinus Bleeding:

FLOSEAL Matrix has been used as a hemostatic agent for the control of operative and post-operative bleeding (epistaxis) during nasal/sinus surgery in 18 patients (30 application sites). Patients were followed for 24 hours following surgery and all complications and episodes of epistaxis were recorded during this period. Intraoperative bleeding stopped in 30 of 30 (100%) application sites. No intraoperative complications were reported in this group. One patient presented with epistaxis 6 hours postoperatively; this patient was treated uneventfully and released from the hospital on the first postoperative day.

Use of FLOSEAL Matrix as a Hemostatic Agent in Cardiovascular Surgery

FLOSEAL Matrix was tested in a prospective randomized controlled trial¹ in 209 patients and compared to 206 control patients treated with either oxidized regenerated cellulose SURGICEL hemostat or purified porcine skin gelatin GELFOAM[®] hemostat. Patients underwent elective cardiac and/or thoracic aortic operations.

Study endpoints included rate of successful intraoperative hemostasis, time required for hemostasis (defined as operative time comprised between decannulation and closure of the sternum), overall postoperative bleeding, rate of transfusion of blood products, rate of surgical revision for bleeding, postoperative morbidity and intensive care unit stay. Patients were followed-up for 96 hours postoperatively.

In the cohort of patients having intraoperative bleeding treated with FLOSEAL Matrix (n=110), a decreased hemostasis time was observed (32.1 ± 5.4 minutes vs. 56.3 ± 8 minutes, p <0.001) as compared to control (n=104). Furthermore a decrease in transfusion of blood products (28.2% vs. 60.6%, p <0.001), revision rate (4.5% vs. 13.5%, p=0.04) and minor complications (20.9% vs. 33.6%, p = 0.04) was observed in the FLOSEAL Matrix group as compared to control. Minor complications were defined as either renal failure, respiratory insufficiency or inotropic support lasting more than 24 hours.

No difference in major complications (stroke, shock, sepsis or myocardial infarction) or ICU stay was observed between groups.

Reference

- Nasso G, Piancone F, Bonifazi R, et al. Prospective, randomized clinical trial of the FloSeal matrix sealant in cardiac surgery. *Ann Thorac Surg* 2009;88(5):1520-1526.

How Supplied:

FLOSEAL Matrix is provided in the configuration shown in the table below.

FLOSEAL Hemostatic Matrix Kit Configuration	
Gelatin Matrix Component	Thrombin Component
<ul style="list-style-type: none"> 1 x 5 mL Syringe with Gelatin Matrix 1 x 5 mL syringe for Matrix Preparation with integral female Luer connector 1 x Bowl for Thrombin 1 x "Thrombin" sticker for Bowl 2 x Applicator tips 	<ul style="list-style-type: none"> 1 x vial Thrombin (Human), Vapor Heated, Solvent/Detergent treated, 2500 units 1 x Pre-filled 0.9% Sodium Chloride Solution syringe, 5mL, USP Injection 1 x Needle-free vial adapter 1 x "Thrombin" sticker for syringe

The package also includes this FLOSEAL Matrix Instructions for Use.

Directions for Use:

Thrombin must be added to the Gelatin Matrix prior to use.

Inspect the integrity of the contents of the FLOSEAL Matrix kit. If the packaging or vials have been damaged or opened, do not use.

Opening the Package

- Place the Thrombin vial, pre-filled Sodium Chloride Solution syringe, and the vial adapter outside of the sterile field.
- Open the outer package containing the Gelatin Matrix Component and deliver the sterile inner package to the sterile field. Once transferred to the sterile field, the inner package may be opened at any time.

Preparing the Thrombin Solution

- Remove the packaging and Luer cap from the pre-filled Sodium Chloride Solution syringe. Remove the plastic flip-off cap from the Thrombin vial. Disinfect the rubber stopper of the Thrombin vial with a germicidal solution and allow to dry. Do not use iodine-containing preparations such as betadine for disinfection.
- Remove the lid from the vial adapter packaging. While gripping the vial adapter packaging, attach the pre-filled Sodium Chloride Solution syringe to the Luer connector of the vial adapter and remove the remaining packaging from the vial adapter.
- While holding the vial adapter, pierce the rubber stopper of the Thrombin vial. Transfer the entire contents of the Sodium Chloride Solution syringe into the Thrombin vial. Leave the syringe attached to the vial adapter and affix the "Thrombin" sticker to it.
- Gently swirl the Thrombin vial with vial adapter and syringe attached until the Thrombin is completely dissolved. Once reconstituted, the Thrombin Solution can be used immediately or may be stored in the vial up to four (4) hours.
- Aspirate the Thrombin Solution into the syringe now labeled "Thrombin".
- The vial may be vented briefly prior to aspiration to facilitate the transfer of the solution in to the syringe. Loosen the syringe connection a quarter (1/4) turn then re-tighten. Ensure that the syringe and vial adapter are firmly connected prior to drawing the solution into the syringe.
- After transfer is complete, disengage the syringe containing Thrombin Solution from the vial adapter.
- Place the Thrombin bowl provided in the Gelatin Matrix Component package onto the sterile field and apply the "Thrombin" sticker to it. Transfer the Thrombin Solution from the syringe into the sterile field by dispensing its contents into the Thrombin bowl.
- Discard the empty Thrombin vial, the vial adapter, and the syringe labeled "Thrombin" appropriately.

Mixing the Thrombin Solution into the Gelatin Matrix

- An empty 5 mL syringe with an integral female Luer connector is provided with the Gelatin Matrix Component. Using this syringe, aspirate the Thrombin Solution from the Thrombin bowl into this empty 5 mL syringe to the indicated mark (4 mL).
- Remove the Luer cap from the Gelatin Matrix Syringe carefully to avoid spilling the Gelatin Matrix granules. Connect the Gelatin Matrix Syringe to the syringe containing the Thrombin Solution. Push the plunger of the Thrombin Solution syringe to fully pass the solution quickly into the Gelatin Matrix Syringe. This constitutes "one pass". Transfer the Gelatin Matrix-Thrombin Solution mixture back and forth between the syringes for a total of at least 20 passes. While starting to mix, do not try to force large, dry clumps of the Gelatin Matrix through the Luer connector, as it may clog. After the first several passes, most of the Gelatin Matrix should be hydrated, and the contents should then be passed rapidly between the syringes to promote thorough mixing. The FLOSEAL Matrix material should be in the syringe labeled FLOSEAL at the completion of mixing.
- Allow 30 seconds after preparation before product is applied to help ensure optimal product consistency and performance. To prevent premature drying of FLOSEAL Matrix, syringes can be kept connected until product is required.
- If desired, connect an Applicator tip to the FLOSEAL Matrix syringe. FLOSEAL Matrix may also be extruded directly from the syringe.
- FLOSEAL Matrix may be used up to eight (8) hours after mixing with the Thrombin Solution.

FLOSEAL Matrix Placement/Application Steps:

Do not inject FLOSEAL Matrix into blood vessels. See the "Contraindications," "Warnings," "Precautions," and "Adverse Events" sections contained in these Instructions for Use.

For best results, FLOSEAL Matrix should be in complete contact with the actively bleeding tissue surface.

The particles of FLOSEAL Matrix swell approximately 10-20% upon contact with blood or other fluids. Maximum swell volume is achieved within about 10 minutes.

Application Technique

Key Application Steps

- Apply FLOSEAL Matrix directly to the source of bleeding.
- Maintain FLOSEAL Matrix at the site (source of bleeding) for two (2) minutes with gentle approximation.
 - For ease of application, steps 1 and 2 may be combined by applying FLOSEAL Matrix onto a moist gauze or cottonoid and use it to deliver FLOSEAL Matrix directly to that source of bleeding.
- Use adequate amounts of FLOSEAL Matrix to completely cover the tissue defect.
- Work quickly.
- Always irrigate excess FLOSEAL Matrix away gently so as not to disturb the new clot and use suction to remove excess.







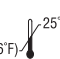


Detailed Application Steps

- Identify the source of bleeding at the tissue surface. This is the target site for FLOSEAL Matrix application.
- Manually approximate a gauze sponge moistened with sterile (non-heparinized) saline against the bleeding surface and use the Applicator tip (or syringe tip) to dispense FLOSEAL Matrix between the sponge and the bleeding surface. The gauze sponge will hold FLOSEAL Matrix in place against the bleeding surface in the presence of active bleeding. Apply enough FLOSEAL Matrix to create a small "mound" of material at the source of bleeding.
- For tissue defects ("divots" or "craters"), begin applying FLOSEAL Matrix at the deepest part of the lesion, and continue applying material as the syringe (or Applicator tip, if used) is withdrawn from the lesion. This "back-filling" action will ensure that FLOSEAL Matrix comes into contact with the entire bleeding surface at the tissue defect.
- Apply a gauze sponge moistened with sterile (non-heparinized) saline to approximate the FLOSEAL Matrix against the bleeding surface, conforming it to the lesion.
- After approximately two (2) minutes, lift the gauze sponge and inspect the wound site. Once bleeding has ceased, excess FLOSEAL Matrix (not incorporated in the hemostatic clot) should always be removed by gentle irrigation and suctioned away from the treatment site.
- If the gauze sponge adheres to the newly formed clot, irrigate with non-heparinized saline to minimize disruption of the clot.
- In cases of persistent bleeding, indicated by saturation and bleeding through the granules, insert the Applicator tip through the center of the mass of previously placed FLOSEAL Matrix to deliver fresh FLOSEAL Matrix as close as possible to the tissue surface. After re-application of FLOSEAL Matrix, resume approximation with a gauze sponge for up to another two (2) minutes, and then inspect the site again. Repeat re-application if necessary.
- Once bleeding has ceased, excess FLOSEAL Matrix, material not incorporated in the hemostatic clot, should always be removed by gentle irrigation and suctioned away from the treatment site.
- Do not disrupt the FLOSEAL Matrix clot by physical manipulation. FLOSEAL Matrix incorporated in the hemostatic clot should be left *in situ*.

Storage Conditions:

The FLOSEAL Matrix kit should be stored at 2 - 25°C (36 - 77°F). **Do not freeze.**

Definition of Symbols:

-  Caution: Consult accompanying documents
-  Consult instructions for use
-  Do not reuse
-  Sterile
-  Do not inject into blood vessels
-  Do not use if package is damaged
-  Temperature Limitations
2°C (36°F) to 25°C (77°F)
-  **Rx ONLY**
Caution: Federal Law (United States) restricts this device to sale by or on the order of a licensed healthcare practitioner
-  **REF**
Reorder/catalog No. 1503350 (Case of 6)

Baxter Healthcare Corporation
21026 Alexander Court
Hayward, CA 94545, USA

Baxter and Floseal are trademarks of Baxter International Inc.

Gelfoam is a registered trademark of Pharmacia & Upjohn Company LLC.

Surgicel is a registered trademark of Johnson & Johnson Corporation.

Patented: Please see www.baxter.com/productpatents/

U.S. - Customer Service 1-888-229-0001

Label Code: 0734989

Rev. Date: 2016-11-01